


InertSearch™ for LC

Inertsil® Applications

Analysis of 5 Nucleobases (Inertsil CX)

Data No. LA740-0489

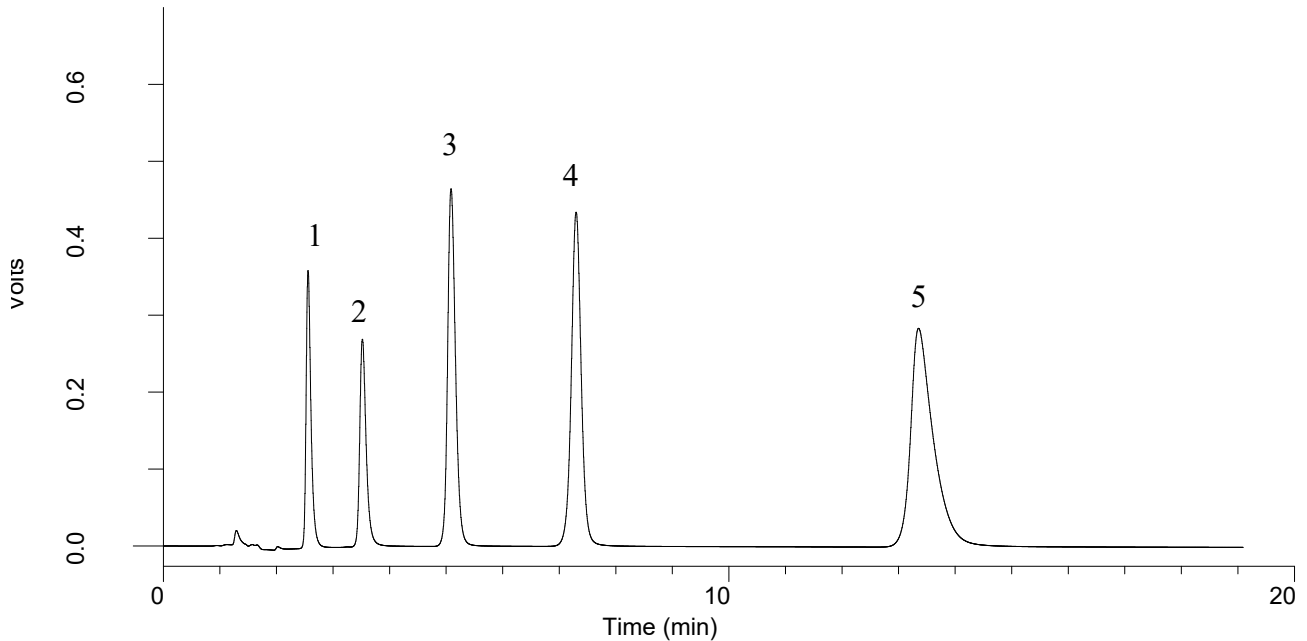

Conditions

Column : Inertsil CX (5 µm, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-07145
Eluent : 0.2 M HCOONH₄ (pH 5.0, HCOOH)
Flow Rate : 1.0 mL/min
Col. Temp. : 40 °C
Detection : UV 260 nm
Injection Vol. : 10 µL

Analyte

1. Uracil (0.05mg/mL)
2. Thymine (0.15mg/mL)
3. Guanine (0.15mg/mL)
4. Cytosine (0.15mg/mL)
5. Adenine (0.15mg/mL)