


InertSearch™ for LC

Inertsil® Applications

Analysis of Vitamin K (Inertsil ODS-2)

Data No. LA765-0696

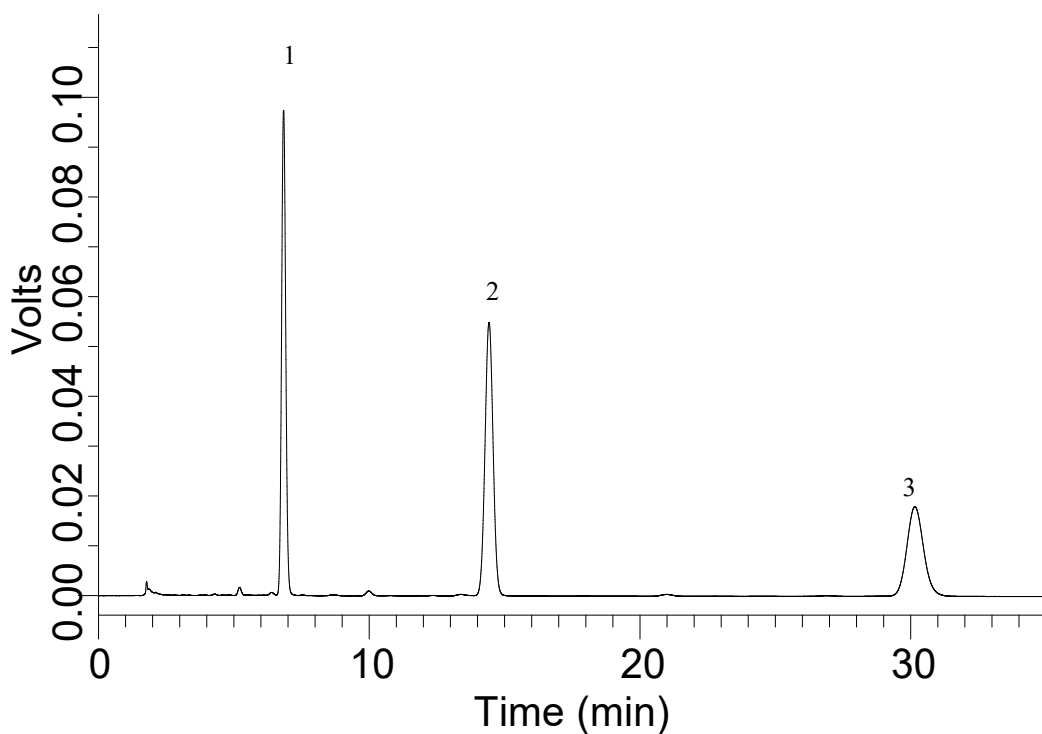

Conditions

System : LC800 HPLC system
Column : Inertsil ODS-2 (5 μ m, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-01101
Eluent : CH₃CN
Flow Rate : 1.0 mL/min
Col. Temp. : 40 °C
Detection : UV 270 nm (LC800 UV Detector)

Analyte

1. Vitamin K2 (MK-4)
2. Vitamin K1
3. Vitamin K2 (MK-7)