


InertSearch™ for LC

Inertsil® Applications

Analysis of Peptides (MonoSpray C18 Nano)

Data No. LA798-0700

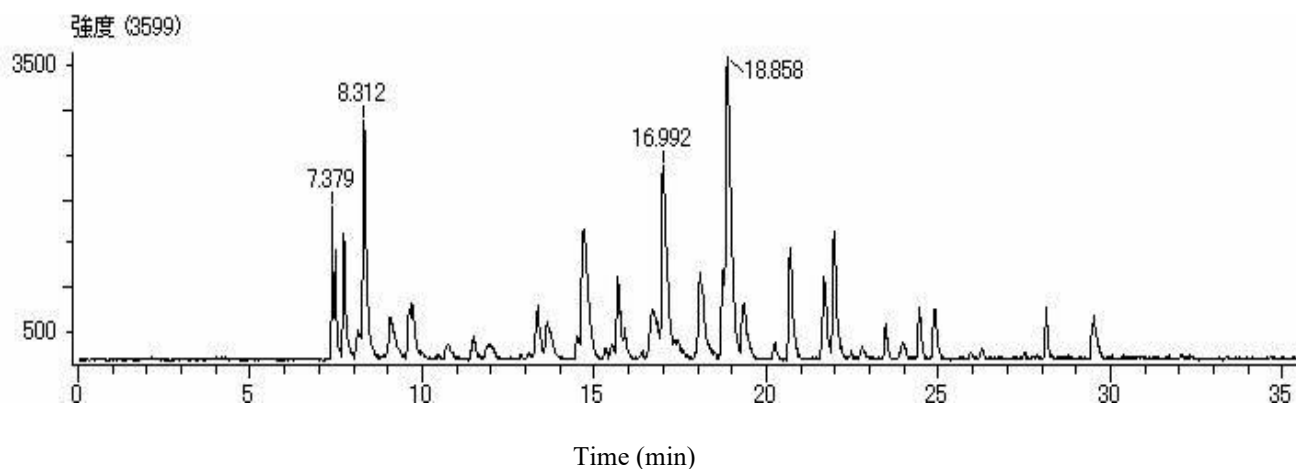

Conditions

System : AccuStream 700 series
Column : MonoSpray C18 Nano (50 x 0.05 mm I.D.)
Column Cat. No. : 5010-20017
Eluent : A) 0.1 % HCOOH in CH₃CN
B) 0.1 % HCOOH
A/B = 10/90 – 30 min – 50/50 – 40 min hold, v/v
Flow Rate : 0.5µL/min
Col. Temp. : Ambient
Detection : Nono-ESI MS positive TIC (JMS-T100-LC, m/z: 500 – 1500)
Injection Vol. : 5 µL (Injected sample was preconcentrated by trap column)

Analyte

Tryptic Digest of BSA