


InertSearch™ for LC

Inertsil® Applications

Analysis of Peptides

Data No. LA837-0700

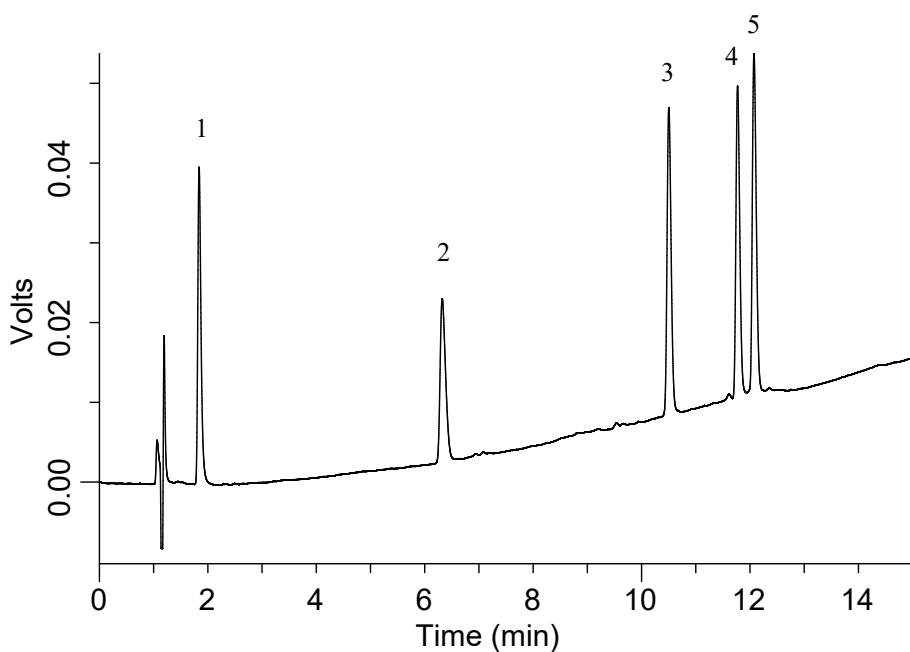

Conditions

System : AccuStream Capillary HPLC system
Column : MonoCap C18 Fast-flow (150 x 0.2 mm I.D.)
Column Cat. No. : 5020-10161
Eluent : A) 0.1 % TFA in CH₃CN
 : B) 0.1 % TFA in H₂O
 : A/B = 10/90 - 15 min - 40/60, v/v
Flow rate : 5 µL/min
Col. Temp. : Ambient
Detection : UV 210 nm (MU701 18 nl cell)
Injection Vol. : 100 nL
Analyte : Peptides

Solutes:

1. Gly-Tyr
2. Val-Tyr-Val
3. Methionine Enkephalin
4. Leucine Enkephalin
5. Angiotensin II