


InertSearch™ for LC

Inertsil® Applications

Analysis of Polycyclic Aromatic Hydrocarbons (Based on the Condition of Japanese Guideline for Particulate Matter 2.5)

Data No. LA933-0871

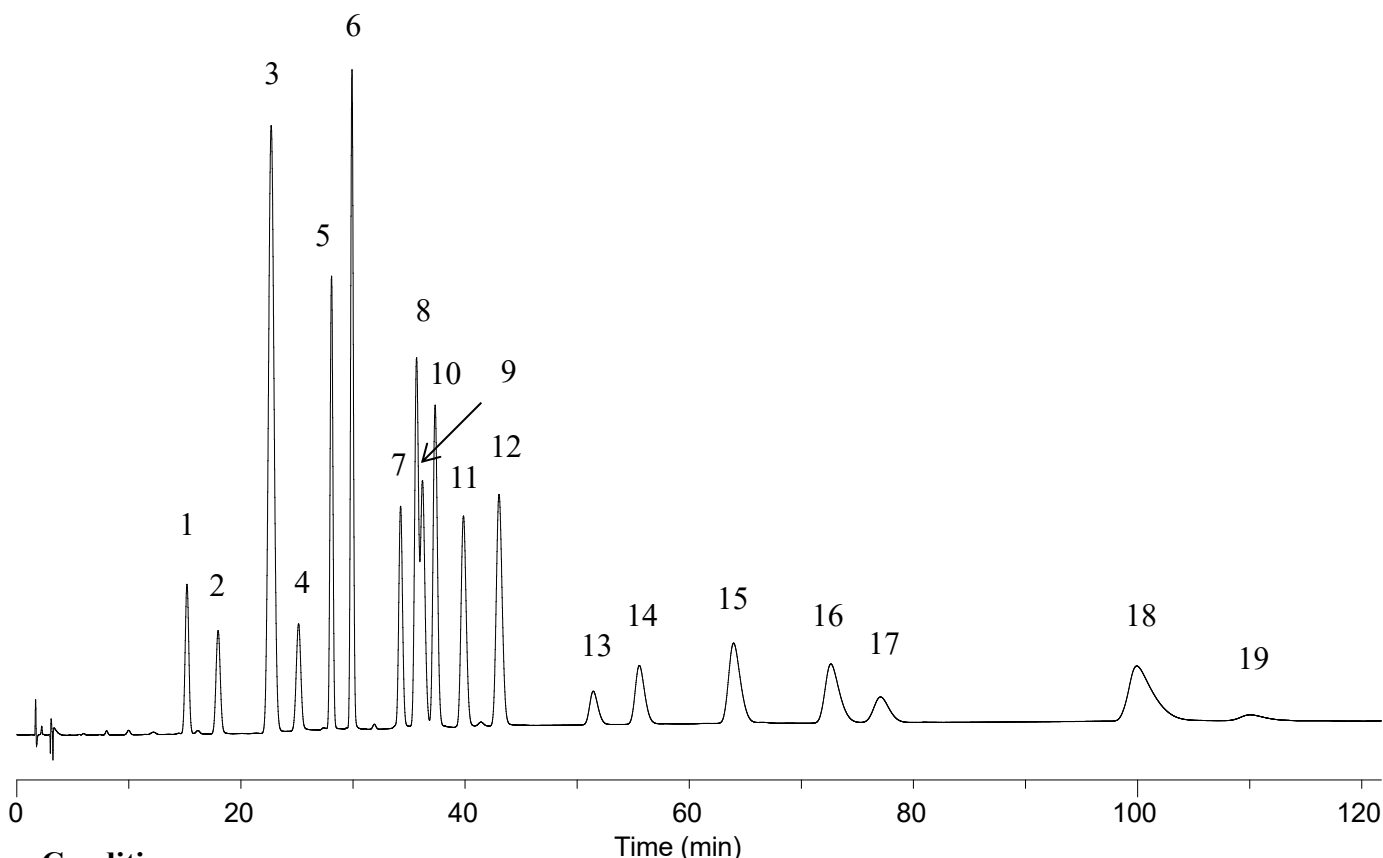

Conditions

System : GL-7400 system
Column : Inertsil ODS-P (5 μm, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-02002
Eluent : A) CH₃CN
 : B) H₂O
 : A/B = 80/20 – 20 min – 80/20 – 5 min
 : – 100/0 – 100 min – 100/0, v/v
Flow rate : 1.0 mL/min
Col. Temp. : 15°C
Detection : UV 254 nm (GL-7452A PDA Detector)
Injection Vol. : 20 μL
Analyte : PAHs

Analyte:

1. Fluoranthene	(10 mg/L)
2. Pyrene	(10 mg/L)
3. Triphenylene	(10 mg/L)
4. <i>p</i> -Terphenyl	(10 mg/L)
5. Benzo[a]anthracene	(10 mg/L)
6. Chrysene	(10 mg/L)
7. Benzo[e]pyrene	(10 mg/L)
8. Benzo[b]fluoranthene	(10 mg/L)
9. Perylene	(10 mg/L)
10. Dibenzo[a,c]anthracene	(10 mg/L)
11. Benzo[k]fluoranthene	(10 mg/L)
12. Benzo[a]pyrene	(10 mg/L)
13. Dibenzo[a,h]anthracene	(10 mg/L)
14. Benzo[g,h,i]perylene	(10 mg/L)
15. Indeno[1,2,3-cd]pyrene	(10 mg/L)
16. Benzo[b]chrysene	(10 mg/L)
17. Dibenzo[a,e]pyrene	(10 mg/L)
18. Picene	(10 mg/L)
19. Coronene	(10 mg/L)