


InertSearch™ for GC

InertCap® Applications

Purity test of Benzyl alcohol

by Japanese Pharmacopoeia

Data No. GA167-0850

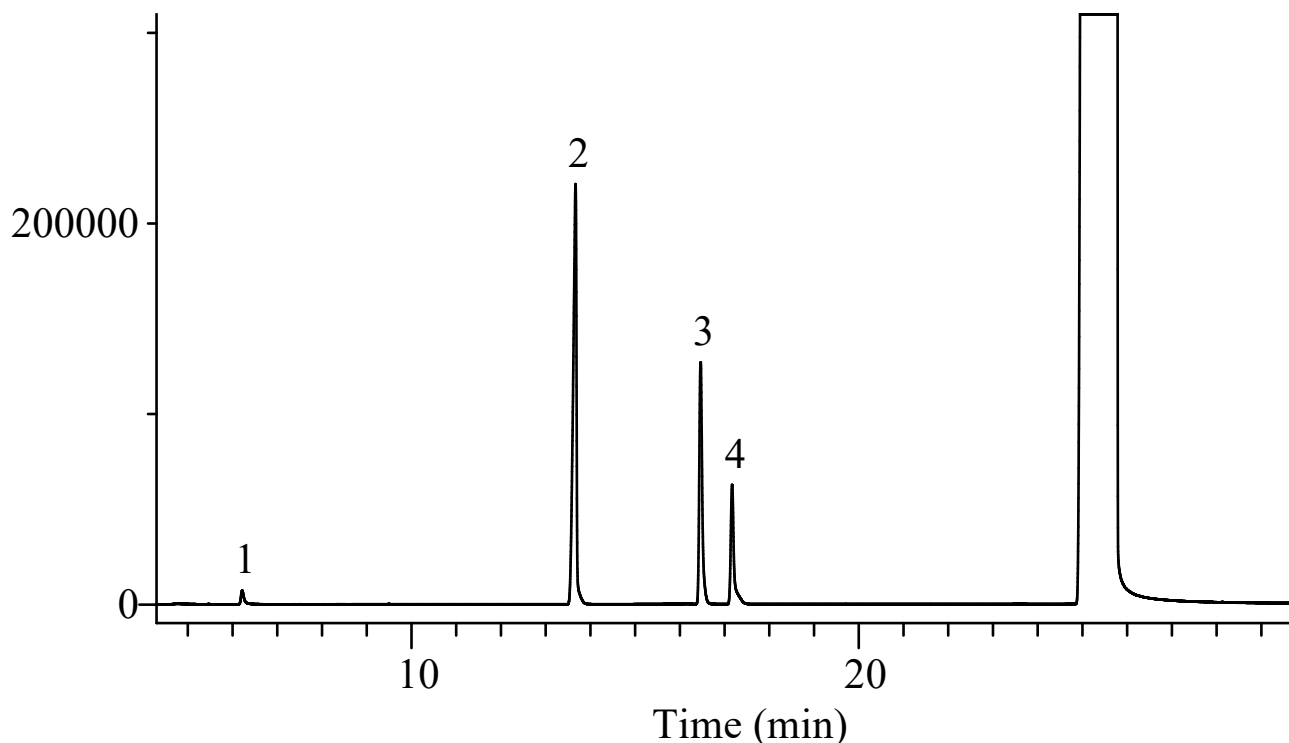

Conditions


System : GC/FID
Column : InertCap WAX
0.32 mm I.D. x 30 m df = 0.50 µm
Col.Cat. No. : 1010-67244
Col. Temp. : 50 °C - 5 °C/min - 220 °C (35 min hold)
Carrier Gas : He 70 kPa
Injection : Splitless 1 min Vent flow 50 mL/min
200 °C
Detection : FID Range 10¹
310 °C
Sample Size : Analyte in Benzyl alcohol
0.1 µL

Analyte :

1. Ethylbenzene 0.1 mg/mL
2. Bicyclohexyl 3.0 mg/mL
3. Benzaldehyde 1.5 mg/mL
4. Cyclohexanemethanol 1.0 mg/mL

*Standard mixture

Chemical Structure

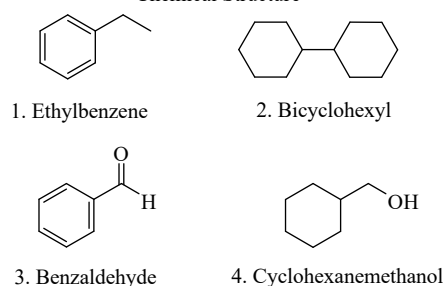

Structures are created using Chemistry 4-D Draw which is provided by ChemInnovation Software, Inc.