


# InertSearch™ for GC

InertCap® Applications

## ε-Caprolactam

Data No. GA171-0904

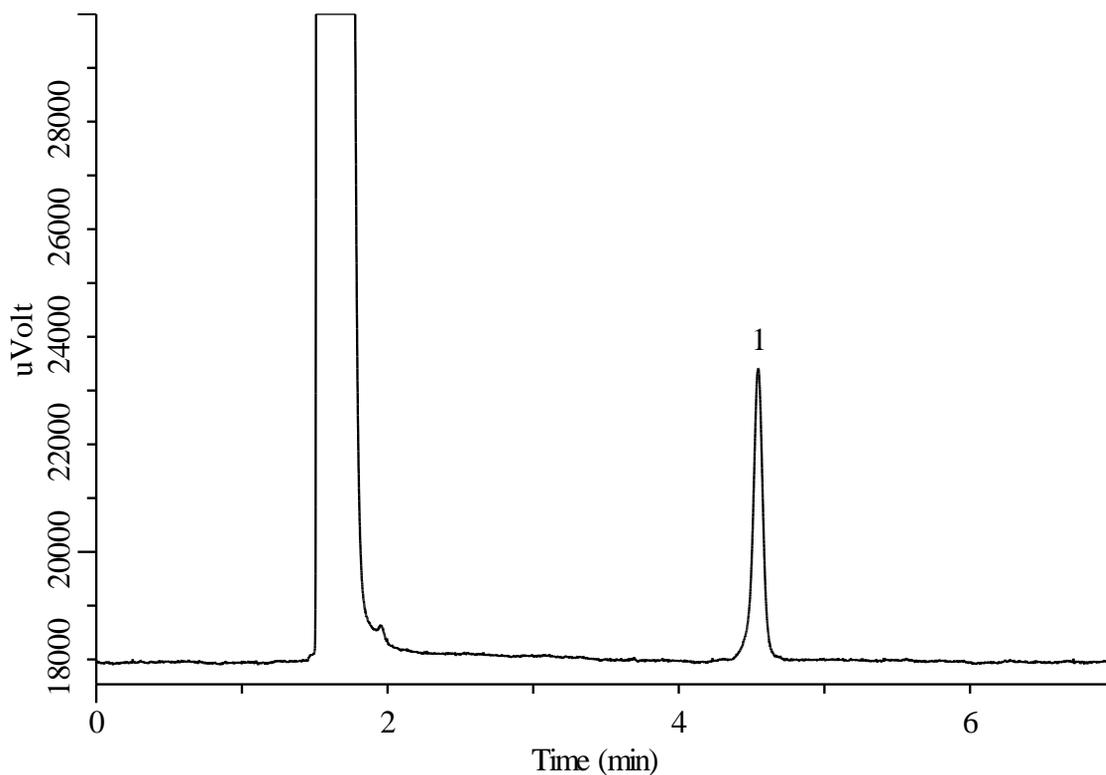

### Conditions

<b>System</b>	: GC/FID
<b>Column</b>	: InertCap 1 0.32 mm I.D. x 30 m df = 5.00 µm
<b>Col.Cat. No.</b>	: 1010-11249
<b>Col. Temp.</b>	: 240 °C
<b>Carrier Gas</b>	: He 110 kPa
<b>Injection</b>	: Split flow 18 mL/min 240 °C
<b>Detection</b>	: FID Range 10 <sup>0</sup> 240 °C
<b>Sample Size</b>	: 15 mg/L in 20% Ethanol solution 1.0 µL
<b>Analyte</b>	: 1. ε-Caprolactam