


InertSearch™ for GC

InertCap® Applications

Hydrocarbons C1-C2

Data No. GA211-0591

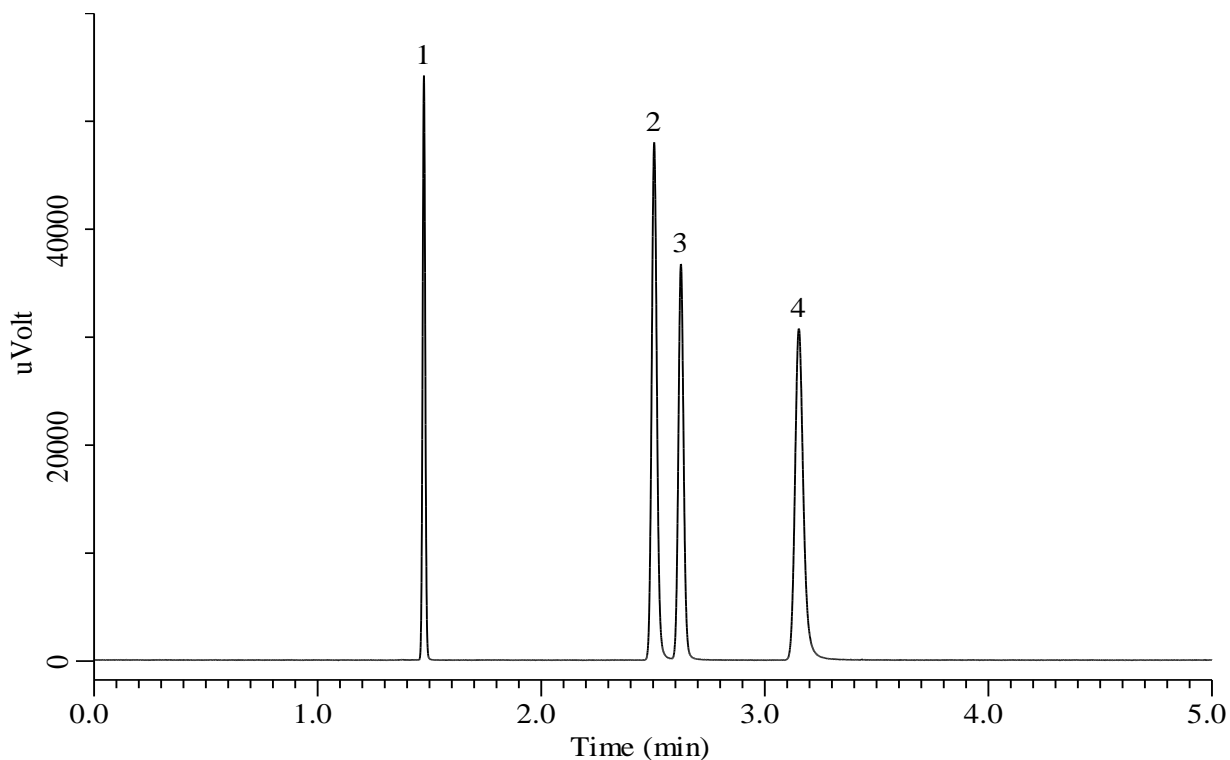

Conditions

System : GC/FID
Column : TC-BOND Q
0.53 mm I.D. x 30 m df = 20 µm
Col. Cat. No. : 1010-40243
Col. Temp. : 40 °C
Carrier Gas : He 30 kPa
Injection : Split flow 50 mL/min
200 °C
Detection : FID Range 10⁰
200 °C
Sample Size : 100 ppm Vol/Vol in He
200 µL

Analyte

1. Methane
2. Ethylene
3. Acetylene
4. Ethane