


InertSearch™ for GC

InertCap® Applications

2-Pyrrolidone

Data No. GA 269-0904

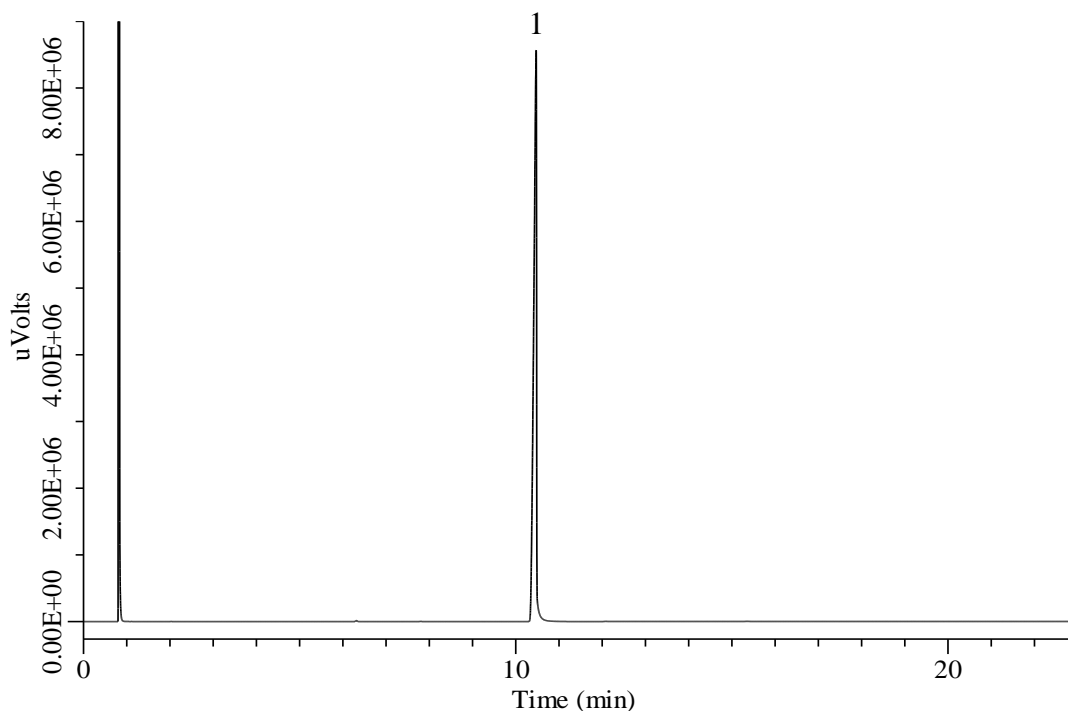

Conditions

System : GC-FID
Column : InertCap Pure-WAX
0.53 mm I.D. x 30 m df = 1.0 µm
Col. Cat. No. : 1010-68445
Col. Temp. : 80 °C(1 min) - 10 °C/min - 190°C(20 min)
Carrier Gas : He 13.0 mL/min Constant Flow
Injection : Split 1:20
200 °C
Injection Vol. : 1 µL
Detection : FID Auto Range
200 °C
Sample : Standard
100 mg/mL in Methanol

Analyte

1. 2-Pyrrolidone