
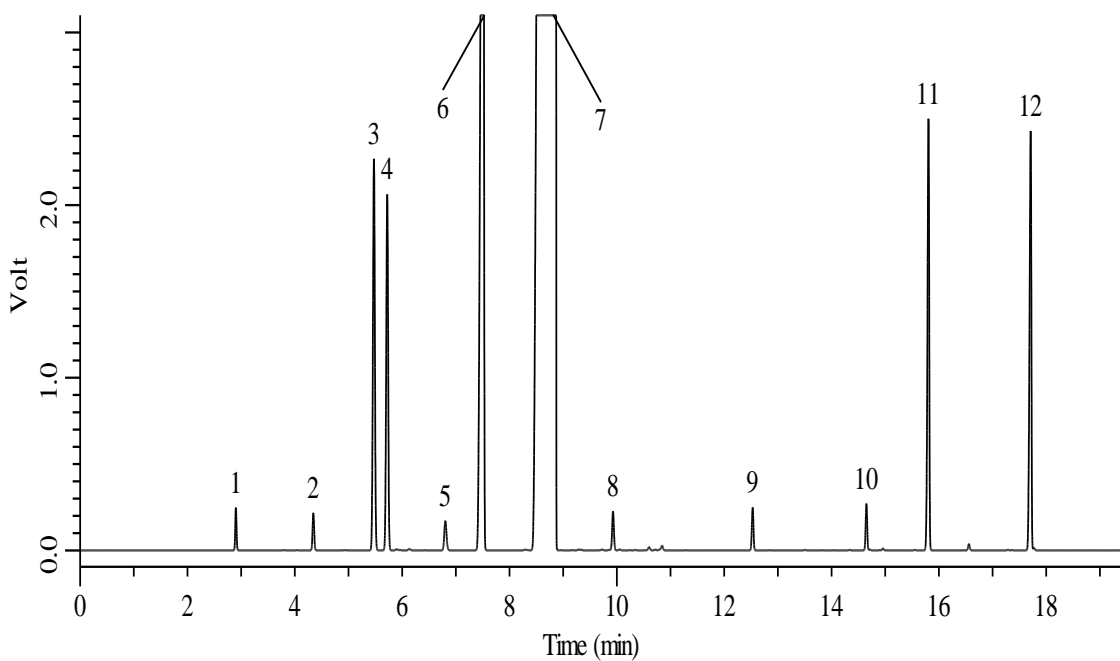


InertSearch™ for GC

InertCap® Applications

Organic hydride

Data No. GA277-0591


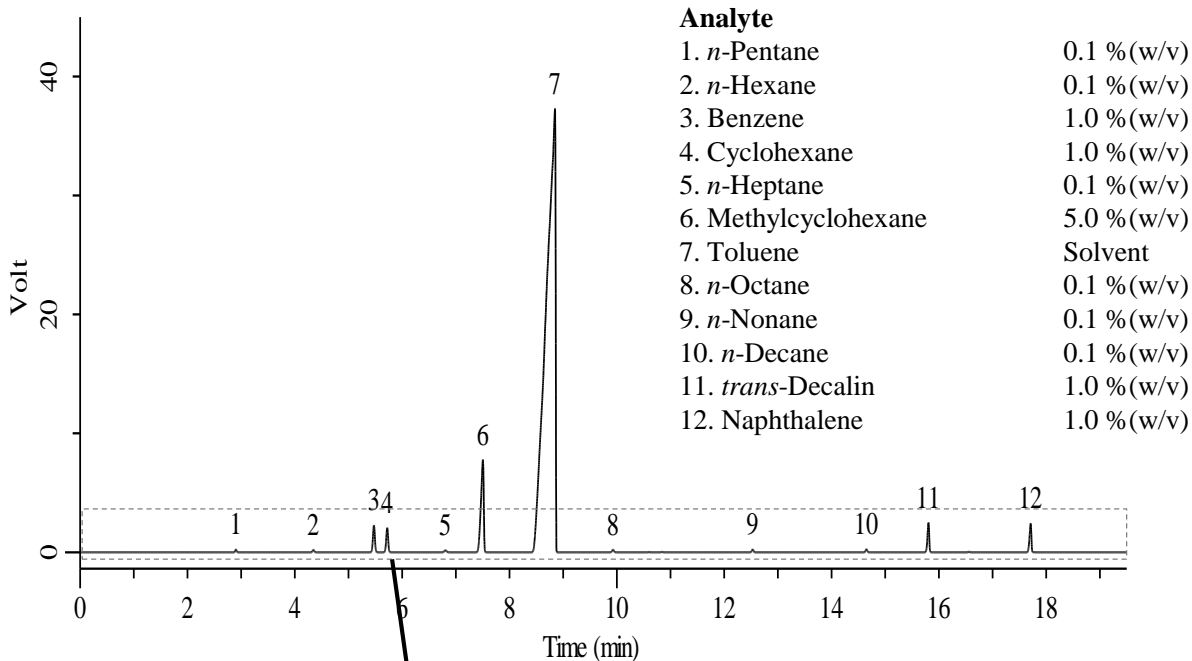


InertSearch™ for GC

InertCap® Applications

Organic hydride

Data No. GA277-0591

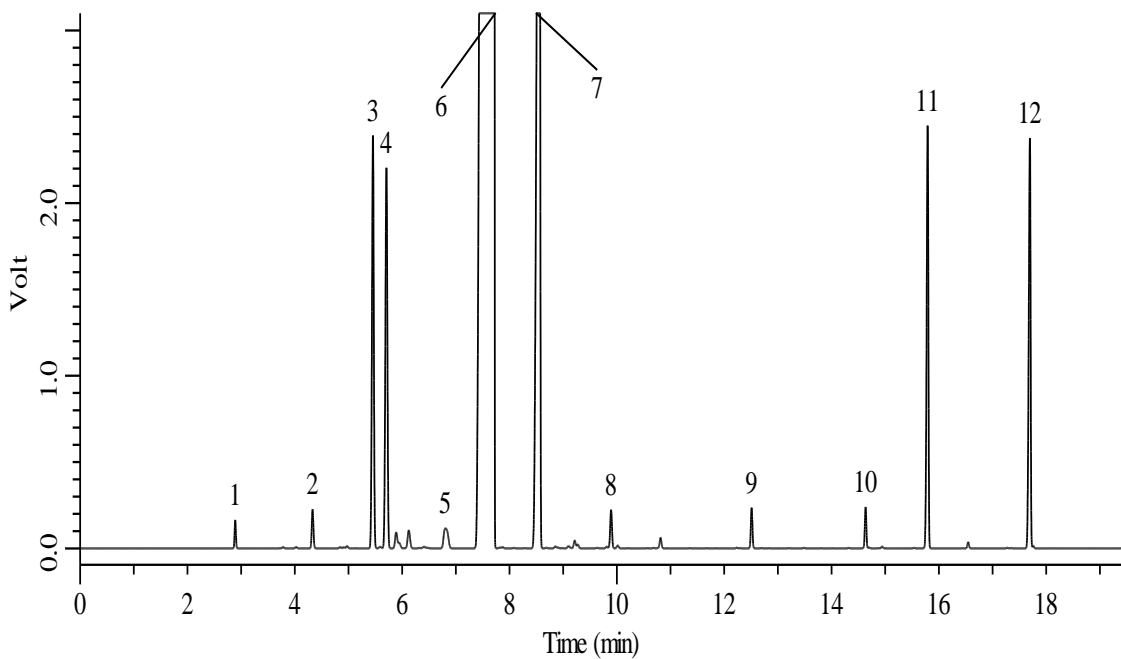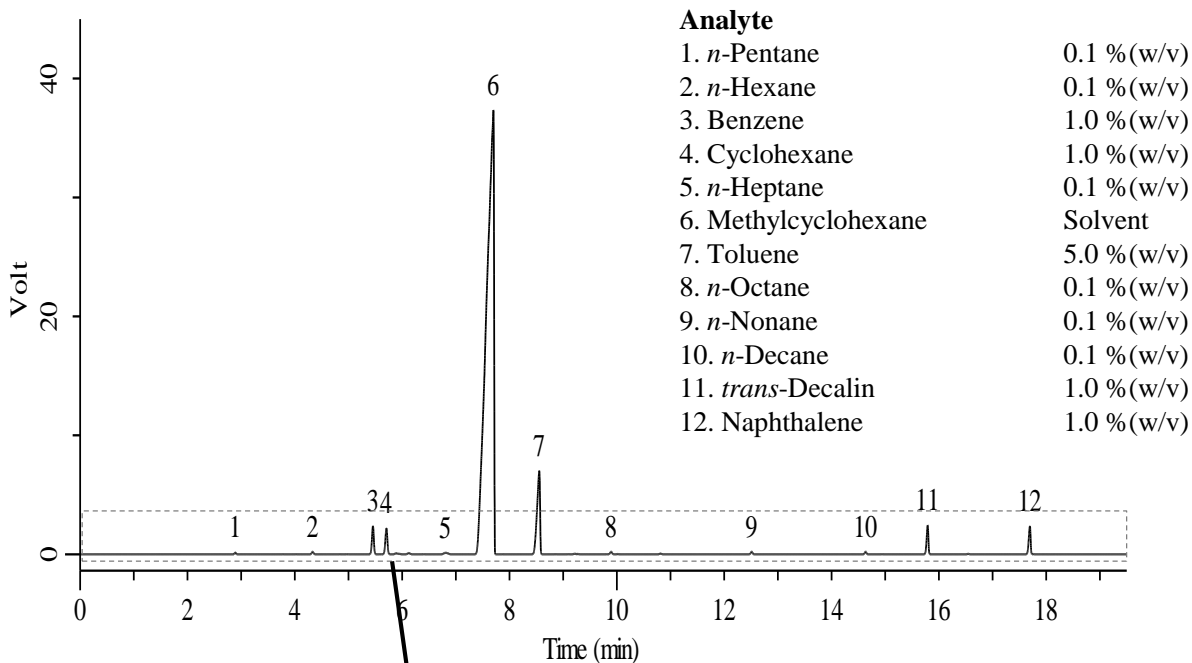

InertSearch™ for GC

InertCap® Applications

Organic hydride

Data No. GA277-0591

Conditions

System : GC-FID
Column : InertCap 1
0.32 mm I.D. x 30 m df = 1.0 µm
Col. Cat. No. : 1010-11245
Col. Temp. : 40 °C (2 min hold) – 5 °C/min – 75 °C – 10 °C/min – 180 °C
Carrier Gas : He 1.7 mL/min
Injection : Split flow 34 mL/min
230 °C
Detection : FID Auto Range
250 °C
Sample Size : 0.5 µL