


Organic Solvents

Data No. GA 284-0644

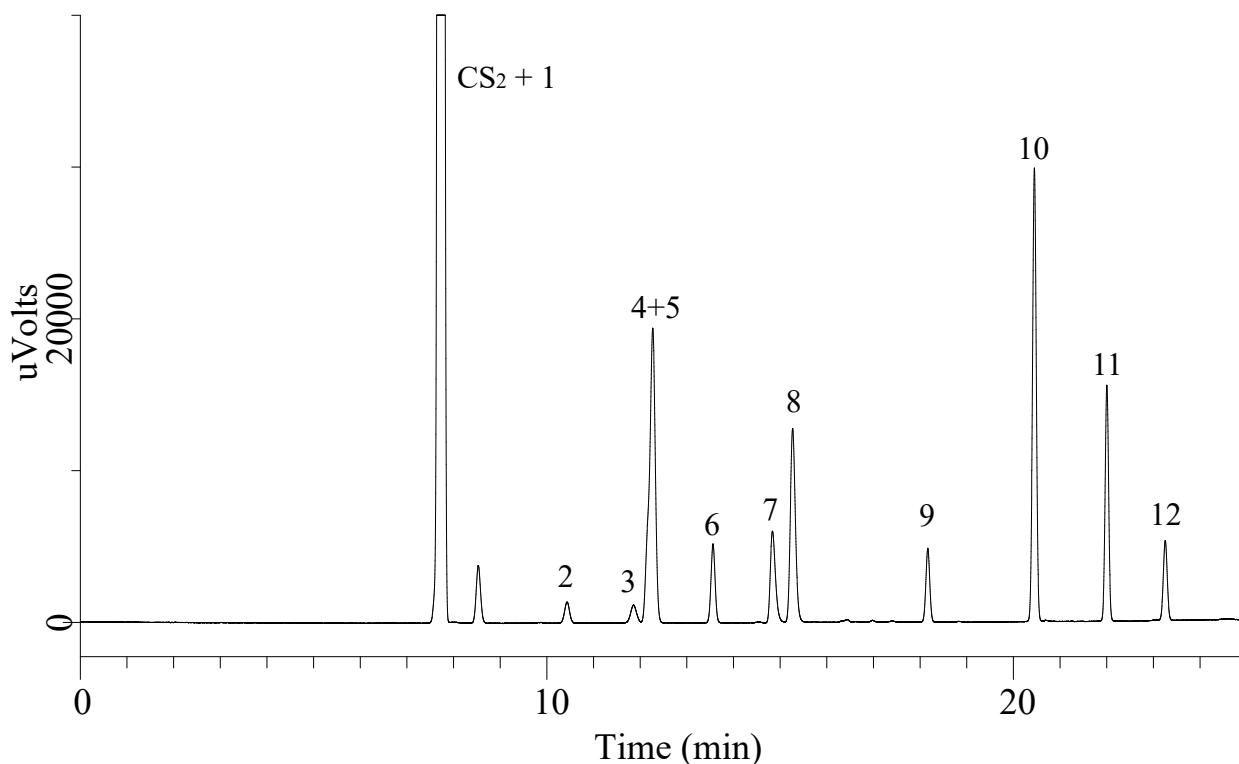

Conditions

System : GC-FID
Column : InertCap AQUATIC
0.53 mm I.D. x 75 m df = 2.0 µm
Col. Cat. No. : 1010-29477
Col. Temp. : 40 °C (5 min hold) - 5 °C/min - 180 °C
Carrier Gas : He 5 mL/min
Injection : Split flow 50 mL/min
220 °C
Detection : FID Auto Range
220 °C
Sample Size : 100 µg/mL in Carbon disulfide
1.0 µL

Analyte : 1. Dichloromethane
2. Chloroform
3. Carbon tetrachloride
4. Benzene
5. 1,2-Dichloroethane
6. Trichloroethylene
7. 1,4-Dioxane
8. 4-Methyl-2-pentanone(MIBK)
9. Tetrachloroethylene
10. Ethylbenzene
11. Styrene
12. 1,1,2,2-Tetrachloroethane