


Organic Solvents

Data No. GA 285-0644

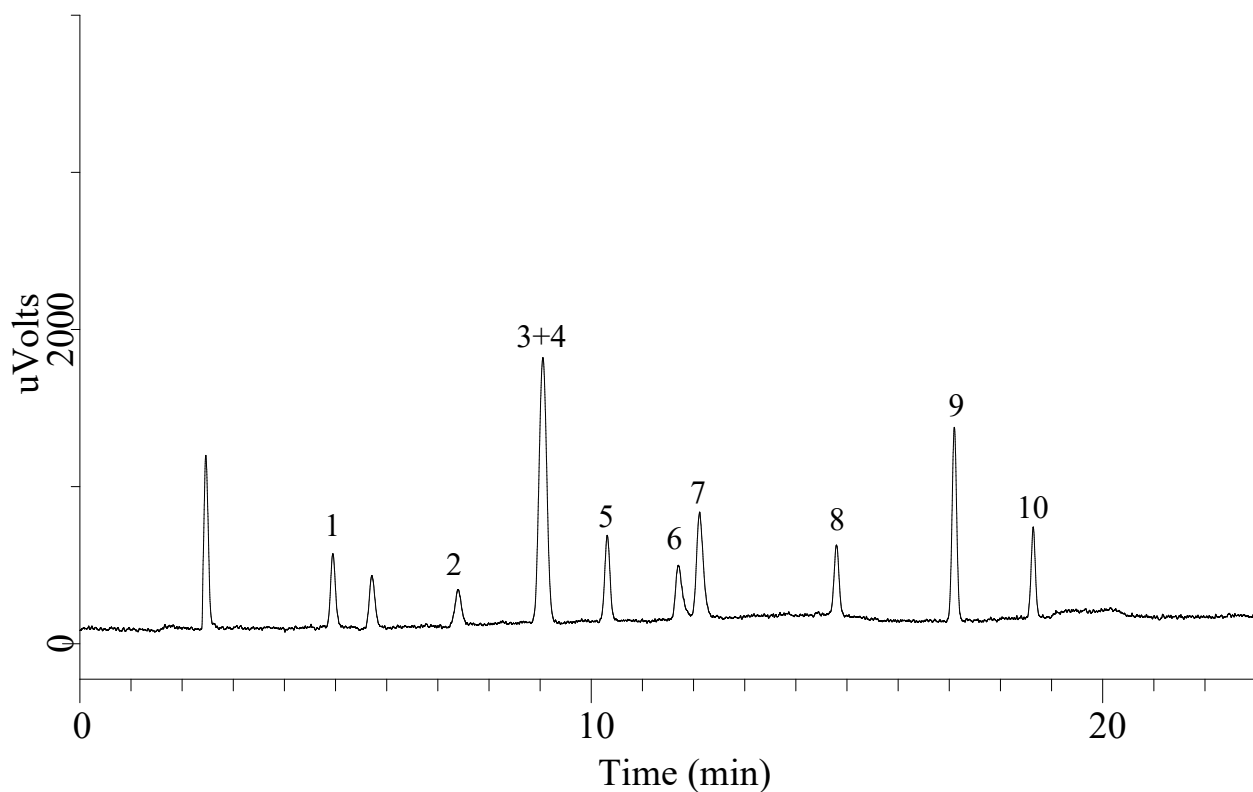

Conditions

System : GC-FID
Column : InertCap AQUATIC
0.53 mm I.D. x 75 m df = 2.0 µm
Col. Cat. No. : 1010-29477
Col. Temp. : 40 °C (5 min hold) - 5 °C/min - 180 °C
Carrier Gas : He 10 mL/min
Injection : Split flow 10 mL/min
220 °C
Detection : FID Auto Range
220 °C
Sample Size : 2-4 ppm(V/V) in N₂
0.5 mL

Analyte : 1. Dichloromethane
2. Chloroform
3. Benzene
4. 1,2-Dichloroethane
5. Trichloroethylene
6. 1,4-Dioxane
7. 4-Methyl-2-pentanone(MIBK)
8. Tetrachloroethylene
9. Ethylbenzene
10. Styrene