


InertSearch™ for LC

Inertsil® Applications

Analysis of Drugs

Data No. LA714-0811

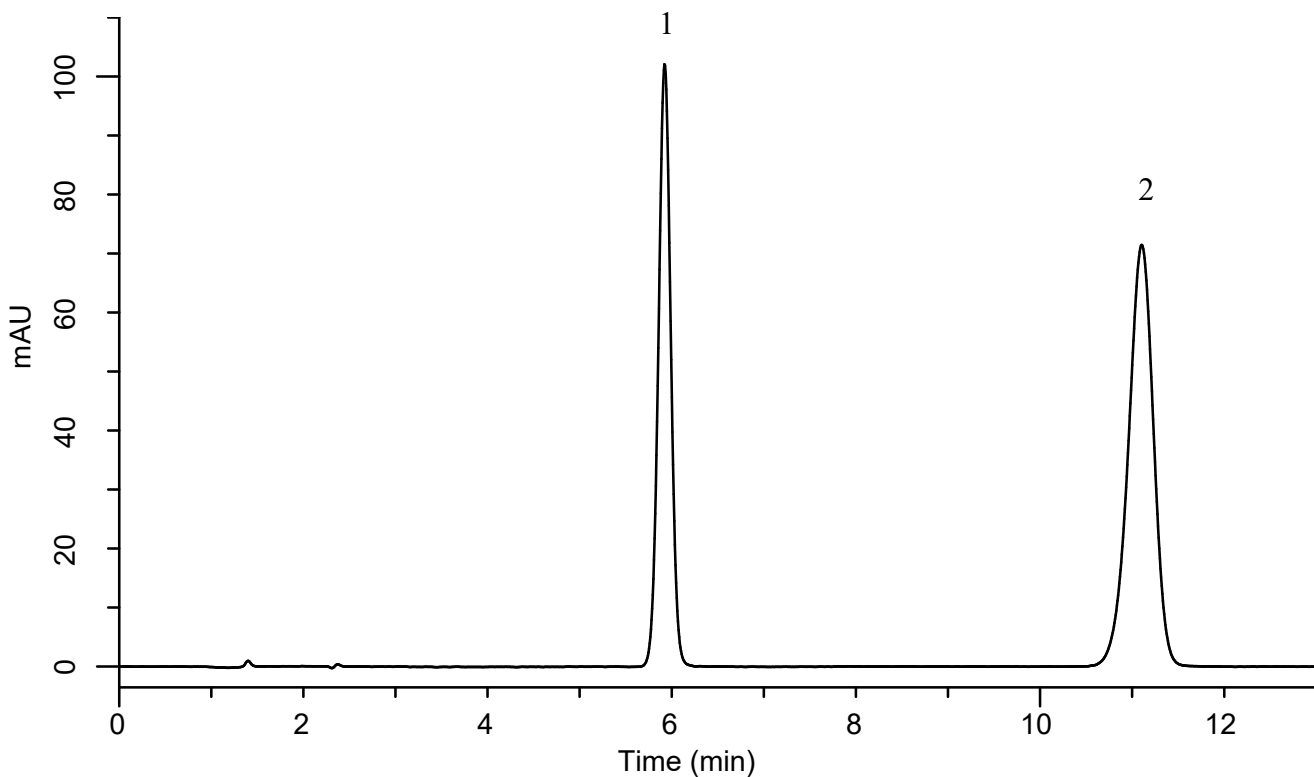

Conditions

System : GL-7400HPLC
Column : Inertsil ODS-3 (5 μ m, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-01731
Eluent : A) CH₃CN
 : B) H₂O
 : A / B = 9 / 91, v / v
Flow Rate : 1.0 mL/min
Col. Temp. : 50°C
Detection : UV 238 nm (GL-7452A PDA Detector)
Injection Vol. : 10 μ L

Analyte

1. Theophylline (100 mg/L)
2. Swertiamarin (100 mg/L)