


# InertSearch™ for LC

Inertsil® Applications

## Analysis of Malachite green and Leucomalachite green

Data No. LA942-0696

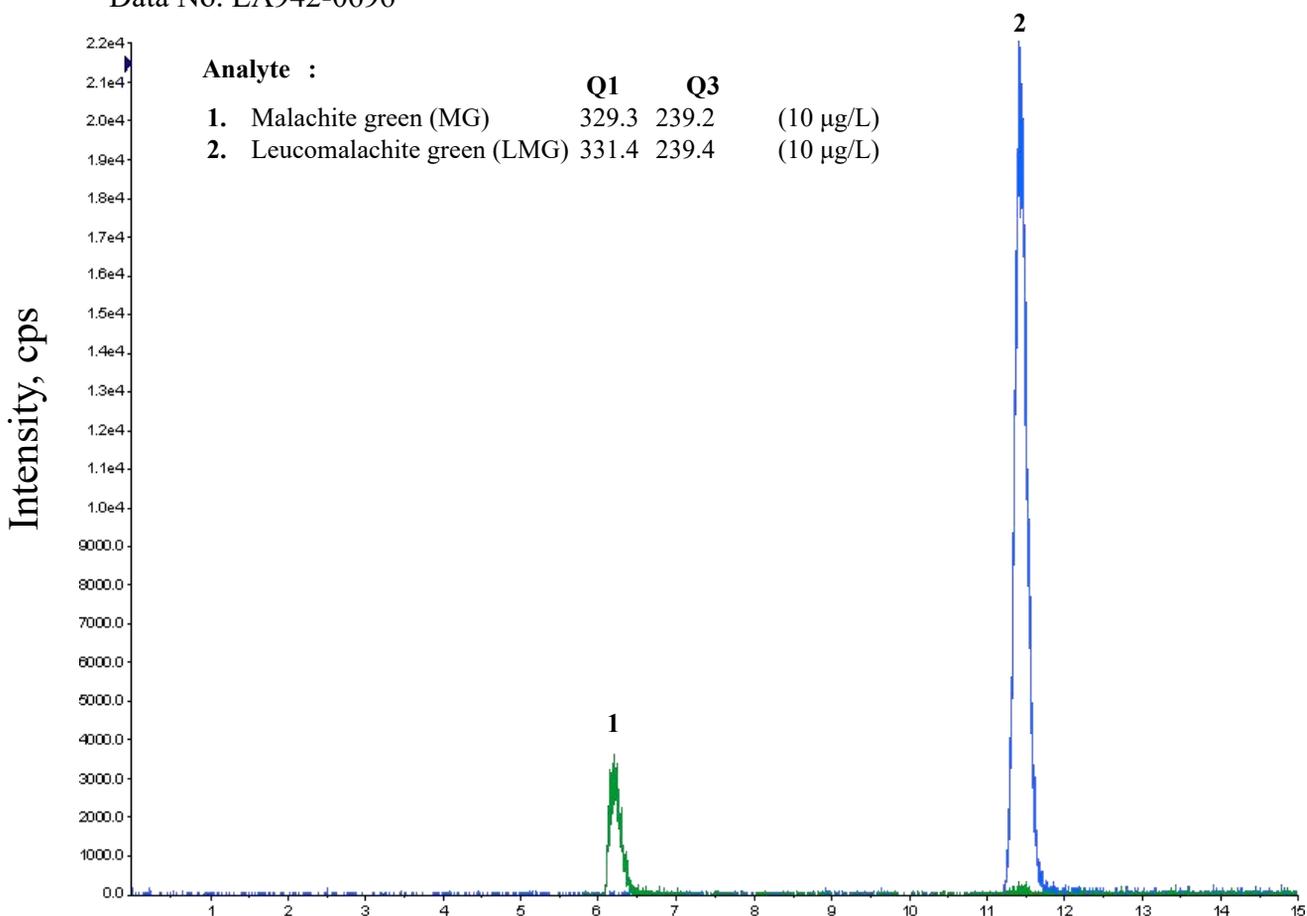

### Conditions

**System** : Symbiosys Pharma system (LC mode)  
API2000®

**Column** : Inertsil ODS-3 (3µm, 150 x 2.1 mm I.D.)

**Column Cat. No.** : 5020-04415

**Eluent** : A) CH<sub>3</sub>CN  
B) 0.1% HCOOH  
A/B = 20/80 -10 min-95/5-10 min -95/5, v/v

**Flow rate** : 0.2 mL/min

**Col. Temp.** : 40 °C

**Detection** : LC/MS/MS (API2000®: ESI, Positive, MRM)

**Injection Vol.** : 10 µL

**Sample** : Standard